
VeriCheck ddPCR™ HEK293 Res DNA Size Kit

User Guide

Catalog #12016813



Table of Contents

Section 1	Introduction	1
Section 2	Kit Contents	1
Section 3	Required Equipment, Reagents, and Consumables	1
Section 4	Sample Preparation: Sample Dilution and DNA Loading	2
Section 5	Reaction Setup	3
Section 6	Thermal Cycling	6
Section 7	Data Acquisition	6
Section 8	Data Analysis	8
	Positive Control–Based Autothresholding	8
	Positive Control–Based Manual Thresholding	9
	Quality Check	9
Section 9	Result Interpretation	11
	Qualitative Assessment	11
	Quantitative Assessment	11
Section 10	Appendices	13
	Appendix A. Performance Specifications	13
	Appendix B. Inhibitors and Extraction-Free Workflow	15
	Appendix C. Understanding Conversion Factor	16
Section 11	Troubleshooting	17
	Too Much Unfragmented DNA	17
	No Negative Droplets	18
	No or Low Total Droplet Counts in Sample Wells Only	18
	No or Low Total Droplet Counts in Both Sample and Control Wells	18
	PCR Inhibition	19
	Mirroring	19
	Droplet Shredding	20
	Incorrect Autothresholding	20

Section 1 Introduction

The VeriCheck ddPCR HEK293 Res DNA Size Kit is designed for the size determination of residual host cell DNA (HCD) using Droplet Digital™ PCR (ddPCR), which provides median size of residual DNA without use of standard curves. HCD carried over during the process of manufacturing therapeutic proteins and vaccines poses safety concerns and must not exceed levels established by regulatory agencies such as the U.S. Food and Drug Administration and the World Health Organization.

Section 2 Kit Contents

This kit uses a three-well test and is sufficient for 100 reactions. Keep at –20°C for long term storage or at 4°C for up to 2 weeks. The kit's contents are listed in Table 1.

Table 1. Kit contents.

Reagent	Vials	Volume
ddPCR Supermix for Residual DNA Quantification (2x)	2	0.75 ml
ddPCR HEK293 Sizing Assay (20x)	1	150 µl
HEK293 Sizing Positive Control	1	150 µl
Nuclease-Free H ₂ O (negative control)	1	1.5 ml

Section 3 Required Equipment, Reagents, and Consumables

Table 2 lists required materials not included in the kit.

Table 2. Additional required materials.

Instrument	Consumables and Reagents*
<ul style="list-style-type: none"> ▪ QX200 Droplet Generator (Bio-Rad™ Laboratories, Inc., catalog #1864002) or Automated Droplet Generator (#1864101) ▪ QX200 Droplet Reader (#1864003) ▪ C1000 Touch Thermal Cycler with 96-Deep Well Reaction Module (#1851197) ▪ PX1 PCR Plate Sealer (#1814000) 	<ul style="list-style-type: none"> ▪ Droplet Generation Oil for Probes (#1863005) or Automated Droplet Generation Oil for Probes (#1864110) ▪ ddPCR Droplet Reader Oil (#1863004) ▪ ddPCR 96-Well Plates (#12001925) ▪ DG8 Cartridges (#1864008) or DG32 Automated Droplet Generator Cartridges (#1864108, 1864109) ▪ ddPCR Buffer Control for Probes (#1863052) ▪ Pipet Tips for AutoDG™ System (#1864120, 1864121) ▪ PCR Plate Heat Seal (#1814040)
<ul style="list-style-type: none"> ▪ QX ONE™ Droplet Digital PCR System (#12006536) ▪ PX1 PCR Plate Sealer (#1814000) 	<ul style="list-style-type: none"> ▪ QX ONE Droplet Generation Oil for Probes (#12006058) ▪ QX ONE Droplet Reader Oil (#12006057) ▪ GCR96 Cartridges (#12006858, 12006859) ▪ ddPCR Buffer Control for Probes (#1863052) ▪ GCR96 Foil Seal (#12006843)

* Adjustable pipets (Rainin or Eppendorf) can be used with the consumables and reagents listed.

Section 4

Sample Preparation: Sample Dilution and DNA Loading

Most samples follow an extraction-free workflow and do not need to be processed using a sample preparation kit. Samples containing high levels of ddPCR inhibitors or background genomic DNA may require dilution prior to setting up the reaction mix. Table 3 contains information regarding suggested dilution ranges for samples that typically occur in different phases of the production process. Examples of early-workflow samples are cell lysate and cell culture media. Examples of mid-workflow samples would be samples in purification stages. Finished products such as chimeric antigen receptor (CAR) T cells would be late-workflow samples. The dilution range of adeno-associated virus (AAV) samples is dependent on the storage buffer and concentration of AAV particles in the starting sample. To determine the optimal dilution for a particular sample, the user can start with the ranges in Table 3 and do a series of dilutions until a unique range for the sample has been determined. See Appendix B for detailed information on allowable concentrations for common inhibitors.

Table 3. Dilution recommendations.

Workflow Stage	Sample Type	Dilution
Early	Cell culture media, cell lysate	3–30x
Middle	Wash buffers used in purification	4–20x
Late	CAR T cells in PBS and HSA	3x
	AAV in PBS	3x
	AAV in water	1x

AAV, adeno-associated virus; CAR, chimeric antigen receptor; HSA, human serum albumin; PBS, phosphate buffered saline.

High levels of cells or unfragmented DNA (including non-HEK DNA) can act as an inhibitor and may require restriction enzyme digestion or dilution. Guidance on DNA loading is listed in Table 4.

Table 4. DNA loading guidance. HEK DNA refers specifically to the HEK DNA to be detected by this kit. Total DNA refers to the total amount of DNA in the sample, including normal human (non-HEK) DNA.

Sample Type	DNA Load	Equivalent Number of Cells
HEK DNA	25 pg to 300 ng	3 to 3e4 cells (triploid cells)
Undigested total DNA	Up to 66 ng	Up to 1e4 cells (diploid cells)
Digested* total DNA	Up to 1 µg digested	Up to 1.7e5 cells (diploid cells)

* Add 0.5 µl HindIII restriction enzyme (at 10,000 units/ml) to each 20 µl ddPCR reaction to digest DNA within droplets. In this case, add only 7.5 µl sample to ensure a total volume of 20 µl per reaction.

Section 5

Reaction Setup

- Determine the plate layout prior to setting up the reaction mix. Each sample and control should have three replicate wells that can be merged to form a complete three-well test. An example plate configuration is shown in Figure 1.

Note: Any unused well in a column from which droplets will be generated must be filled with ddPCR Buffer Control for Probes prior to droplet generation.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Pos Ctrl	Pos Ctrl	Pos Ctrl	#7	#7	#7						
B	Sample #1	Sample #1	Sample #1	#8	#8	#8						
C	Sample #2	Sample #2	Sample #2	#9	#9	#9						
D	Sample #3	Sample #3	Sample #3	#10	#10	#10						
E	Sample #4	Sample #4	Sample #4	#11	#11	#11						
F	Sample #5	Sample #5	Sample #5	#12	#12	#12						
G	Sample #6	Sample #6	Sample #6	#13	#13	#13						
H	Neg Ctrl	Neg Ctrl	Neg Ctrl	#14	#14	#14						

Fig. 1. Example plate configuration for controls and 14 samples. Columns 7–12 can be used to run additional samples. **Three replicate wells** should be included for every sample and control.

- Sample:** DNA sample is added to the ddPCR reaction.
- Pos Ctrl (HEK293 Sizing Positive Control):** Instead of sample, positive control is added to the ddPCR reaction.
- Neg Ctrl (Nuclease-Free H₂O [negative control]):** Instead of sample, nuclease-free water is added to the ddPCR reaction.

- Thaw all kit components to room temperature. Mix all tubes thoroughly by vortexing at **maximum speed for 15 sec** to ensure homogeneity. Centrifuge briefly to collect contents at the bottom of each tube.

Note: Thorough vortexing (maximum speed for 15 sec) of each component at this step is essential.

- Ensure samples are thoroughly mixed. If necessary, prepare samples to the desired concentration before setting up the reaction mix. The recommended input DNA is 25 pg–300 ng per well. See Appendix B for additional sample dilution recommendations based on common ddPCR inhibitors.
- Assemble a master mix containing supermix, sizing assay, and water, according to the master mix setup guidelines in Table 5. Vortex the master mix thoroughly at maximum speed for 15 sec.

Table 5. Preparation of master mix containing supermix, sizing assay, and water.

Component	Volume per Reaction, μ l		Volume per 96 Reactions + Overage*, μ l
	1 Well	1 Well + Overage*	30 Samples, 3 Wells Each**
ddPCR Supermix for Residual DNA Quantification	10	11	1,320
HEK293 Sizing Assay	1	1.1	132
Nuclease-Free H ₂ O (negative control)	1	1.1	132
Total volume	12	13.2	1,584

* Prepared with overage to account for liquid loss during pipetting steps.

** Volumes in this column are for a total of 96 wells (30 samples [3 wells each], positive control [3 wells], and nuclease-free H₂O/negative control [3 wells]).

5. Prepare reaction mix for each set of three replicate wells by combining master mix (from step 4) and the sample or control.
 - a. Combine master mix and sample or control according to the volumes given in Table 6 for each reaction. Mix each tube thoroughly by vortexing at maximum speed for 15 sec, then centrifuge briefly to collect the contents at the bottom of each tube.

Tip: If several reaction mixes are being prepared at once, they can be made in a column of a ddPCR 96-Well Plate (mixed thoroughly) and the reaction mix from each well of the column can be transferred simultaneously to the neighboring columns using a multichannel pipet.

Table 6. Preparation of reaction mix for three replicate wells.

Component	Volume per Reaction, μ l	Volume per Reaction, μ l	Volume per 3 Reactions + Overage*, μ l
	1 Well	1 Well + Overage*	3 Wells + Overage*, μ l
Master mix	12	13.2	48
Sample/positive control/nuclease-free H ₂ O (negative control)	8	8.8	32
Total volume	20	22	80

* Prepared with overage to account for liquid loss during pipetting steps.

6. Transfer the reaction mix or buffer control to the appropriate droplet generation cartridge and generate droplets based on the system in use:
 - a. **For the QX200 Droplet Digital PCR System:**
 - i. Load 20 μ l of each reaction mix into the sample wells of a DG8 Cartridge preloaded in a DG8 Cartridge holder. The orientation of the manual DG8 Cartridge should mirror the way the droplets will be transferred to the final 96-well plate.
 - ii. Add 20 μ l of buffer control to any unused wells in a column from which droplets will be generated.
 - iii. Load 70 μ l of Droplet Generation Oil for Probes into the oil wells. For detailed instructions, refer to the QX200 Droplet Generator Instruction Manual ([10031907](#)).
 - b. **For the QX200 AutoDG Droplet Digital PCR System:**
 - i. Load 22 μ l of each reaction mix into the wells of a ddPCR 96-Well Plate.
 - ii. Add 22 μ l of buffer control to any unused wells in a column from which droplets will be generated. If an entire column is unused, no buffer control is required in those wells.
 - iii. Seal the plate using the PX1 PCR Plate Sealer at 180°C for **5 sec**, allow to cool briefly, vortex at maximum speed for 15 sec, and centrifuge at 1,000 rcf for 1 min. Visually verify that all the liquid is at the bottom of the well.
 - iv. Place the sealed plate in the Automated Droplet Generator and follow instructions in the Automated Droplet Generator Instruction Manual ([10043138](#)).
 - c. **For the QX ONE ddPCR System:**
 - i. Load 20 μ l of each reaction mix into the wells of a GCR96 Cartridge.
 - ii. Droplet generation occurs in sets of two columns (for example, droplets for columns 1 and 2 are generated simultaneously). Add 20 μ l of buffer control to any unused wells in each set of columns from which droplets will be generated.

- iii. Seal the plate using the PX1 PCR Plate Sealer at 180°C for **0.5 sec**, flip the plate 180° and repeat sealing, allow to cool briefly, and centrifuge at 1,150 rcf for **30 sec**.
- iv. Follow subsequent instructions for loading the plate, as specified in the QX ONE Droplet Digital PCR System and QX ONE Software User Guide (10000116512). Use the appropriate thermal cycling conditions specified in Table 7.

Important: When setting up the plate template, first select the wells to include/exclude under the Exclude tab. Then move to the Edit tab and designate the following:

- **Experiment Type:** Residual DNA Quantification (RDQ)
- **Sample Descriptions:** Determined by user
- **Sample Type:** Determined by user (Pos Ctrl/Unknown/NTC)
- **Conversion Factor:** 0.55 (HEK293, HEK293T)
 - See Appendix C for more information
- **Supermix:** ddPCR Supermix for Residual DNA Quantification
- **Assay Type:** Single Target per Channel
- **Target Name(s):** Determined by user
- **Target Type:** Unkn
- **Signal Ch1:** FAM
- **Signal Ch2:** HEX

Click **Apply** and then **Save** to save the template. Click **Start Run**. See Figure 2 for an example.

Section 6

Thermal Cycling

Follow the instructions for thermal cycling based on the system in use.

- **For the QX200 ddPCR System**, after droplet generation with the QX200 Droplet Generator, carefully transfer each column of the droplet emulsions into a clean ddPCR 96-Well Plate using a P50 multichannel pipet. Seal the plate using the PX1 PCR Plate Sealer at 180°C for 5 sec. Proceed to thermal cycling. Use the appropriate thermal cycling conditions specified in Table 7
- **For the QX200 AutoDG ddPCR System**, remove the droplet plate containing ddPCR droplets from the Automated Droplet Generator. Seal the plate using the PX1 PCR Plate Sealer at 180°C for 5 sec. Proceed to thermal cycling. Use the appropriate thermal cycling conditions specified in Table 7
- **For the QX ONE ddPCR System**, thermal cycling is integrated into and sequentially performed by the system itself. Hence, no additional equipment or sample handling is required for this step

Table 7. Cycling conditions.*

Cycling Step		Temperature, °C	Time	Number of Cycles	Ramp Rate
Hold (QX ONE ddPCR System only)		25	3 min	1	2°C/sec
Enzyme activation		95	10 min	1	
Denaturation		94	30 sec	40	
Annealing/extension		60	1 min		
Enzyme deactivation		98	10 min	1	
Hold	QX200 ddPCR System (optional)	4	∞	1	
	QX ONE ddPCR System (required)	25	1 min	1	

* For the C1000 Touch Thermal Cycler, use a heated lid set to 105°C and set the sample volume to 40 µl.

Section 7

Data Acquisition

Follow the instructions for data acquisition based on the system in use.

- **For the QX200 ddPCR System and the QX200 AutoDG ddPCR System, using QX Manager Standard or Regulatory Edition version 1.2 or higher**, place the sealed 96-well plate in the QX200 Droplet Reader and refer to the QX200 Droplet Reader and QX Manager Software Regulatory Edition User Guide ([10000107224](#)) or the QX200 Droplet Reader and QX Manager Software Standard Edition User Guide ([10000107223](#)) for plate setup. See Figure 2 for an example.

Important: When setting up the plate template, first select the wells to include/exclude under the Exclude tab. Then move to the Edit tab and designate the following:

- **Experiment Type:** Residual DNA Quantification (RDQ)
- **Sample Descriptions:** Determined by user
- **Sample Type:** Determined by user (Pos Ctrl/Unknown/NTC)
- **Conversion Factor:** 0.55 (HEK293, HEK293T)
 - See Appendix C for more information
- **Supermix:** ddPCR Supermix for Residual DNA Quantification

- **Assay Type:** Single Target per Channel
- **Target Name(s):** Determined by user
- **Target Type:** Unkn
- **Signal Ch1:** FAM
- **Signal Ch2:** HEX

Click **Apply** and then **Save** to save the template. Click **Start Run**. See Figure 2 for an example.

- **For the QX ONE ddPCR System using QX ONE Software version 1.1 or higher**, data acquisition is integrated into and sequentially performed by the system itself. Hence, no additional equipment or sample handling is required for this step.

Experiment Type
Residual DNA Quantification (RDQ) Apply

Sample Description1: Pos Ctrl Sample Description2: 1 Apply

Sample Description3: Sample Description4:

Sample Type: Pos Ctrl (highlighted) Conversion Factor: 0.550

Supermix: ddPCR Supermix for Residual DNA Quantification

Assay Type: Single Target per Channel Apply

Target Name	Target Type	Signal Ch1	Signal Ch2
HEK293	Unkn	FAM	None
HEK293_200	Unkn	None	HEX

Well Notes: Plate Notes: Apply

Fig. 2. Example of plate setup in the Plate Editor screen. Positive control wells should have the sample type **Pos Ctrl** (yellow rectangle). Remove additional (more than two) targets and channels.

Section 8

Data Analysis

Positive Control–Based Autothresholding

Positive control–based autothresholding generates thresholds for each channel based on the wells marked as having the sample type **Pos Ctrl**. It then applies those thresholds to all currently selected wells. Refer to the QX200 Droplet Reader and QX Manager Software Regulatory Edition User Guide (10000107224) or the QX200 Droplet Reader and QX Manager Software Standard Edition User Guide (10000107223) for detailed instructions about data analysis. QX ONE Software 1.1 is not compatible with positive control–based autothresholding. See Positive Control–Based Manual Thresholding section below.

To threshold a plate with one assay:

1. Ensure that only the positive control wells for the assay are marked as sample type **Pos Ctrl**. If necessary, navigate to the Plate Editor tab and select the positive control wells. Under Sample Type, use the dropdown menu to select **Pos Ctrl** and click **Apply**. See Figure 2 for details.
2. Navigate to the 2D Amplitude tab. Select the wells of interest on the plate, including the control wells.
3. Click the **Auto with Tilt** button, select **Positive Control Wells** at the bottom of the dropdown menu, and click **OK**. This will autothreshold all selected wells based on the designated positive control wells. See Figure 3 for details.

Important: At least one of the selected wells must have sample type **Pos Ctrl** in order to access this method of thresholding.

To threshold a plate with multiple assays:

If multiple assays are run on the same plate, the wells for each assay must be thresholded separately. Mark as **Pos Ctrl** the positive control wells from only one assay at a time. Positive control wells from all other assays on the plate should be marked as sample type **Unknown**. Use the Sample Type dropdown menu and click **Apply** to mark selected wells.

After thresholding all wells in the first assay, return to the Plate Editor tab and change the positive control wells for the first assay to Sample Type **Unknown**. Then repeat the thresholding tips above for the next assay. Repeat this process until all assays are thresholded.

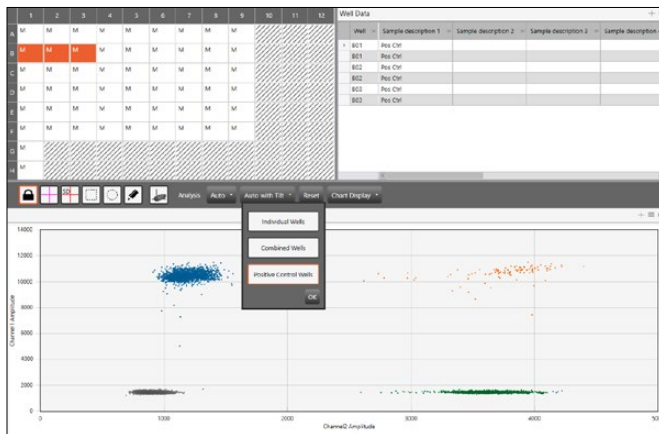


Fig. 3. 2D Amplitude screen when using positive control–based autothresholding with tilt.

Positive Control–Based Manual Thresholding

QX ONE Software version 1.1 does not support autothresholding. To threshold manually, follow the steps below.

1. Navigate to the 2D Amplitude tab and select all sample and control wells of interest on the plate. Click the **Auto with Tilt** button, select **Combined Wells** at the bottom of the dropdown menu, and click **OK**. This will tilt-correct all selected wells and apply a threshold.
2. To manually adjust the thresholds, select the positive and negative control wells for one assay. Select the pink manual thresholding tool on the 2D Amplitude tab and drag the thresholds to adjust them manually.
3. Copy the manual threshold applied to the control wells and paste it into all other wells with the same assay (the sample wells).
4. If multiple assays are run on the same plate, repeat steps 1–3 for each assay.

Tip: Refer to the QX ONE Software Regulatory Edition Instruction Manual ([10000116656](#)) and the Droplet Digital PCR Applications Guide ([bulletin 6407](#)) for detailed instructions and best practices for setting manual thresholds.

Quality Check

1. Navigate to the Event Counts tab.

Recommended best practice: If the droplet count is <10,000 in any well, do not include that well in the final analysis. If all wells fail for a sample, refer to the troubleshooting section.

2. Navigate to the 2D Amplitude tab and examine the control wells for uniformity. Do not use poor-quality positive control well(s) as thresholding controls. See Figures 4 and 5 for example positive and negative control well 2D plots.

Tip: If any poor-quality wells were used to threshold the plate, return to the Plate Editor tab and mark poor-quality positive control well(s) as having the sample type **Unknown**. Then, threshold the plate again using positive control–based autothresholding. Only one positive control well is required for accurate thresholding.

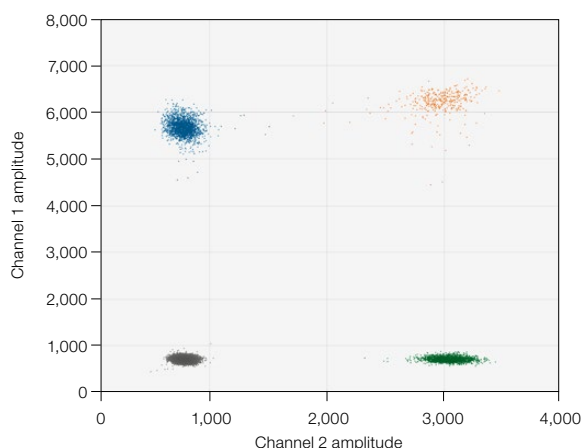


Fig. 4. Example positive control well.

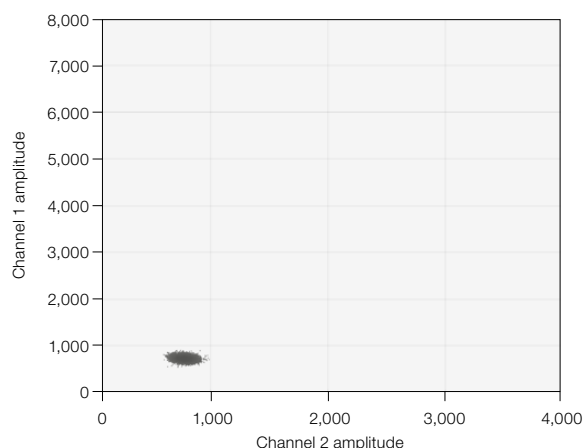


Fig. 5. Example negative control well.

3. Examine the replicate sample wells for uniformity and examine any outliers for problems. Do not include problem wells in the final analysis. If any clusters are cut off by the thresholds (for example, see Figure 16), the thresholds should be manually adjusted. See the troubleshooting section for more details.
4. Navigate to the Concentration tab and examine the control wells and each sample's replicate wells for uniformity. Positive control wells are expected to have a FAM (channel 1) and HEX (channel 2) concentration of 69–275 copies/ μ l (cp/ μ l). Negative control wells are expected to have a FAM and HEX concentration of 0 cp/ μ l.
5. If all replicate wells have similar concentrations, merge the wells for the final analysis. Refer to the QX200 Droplet Reader and QX Manager Software Regulatory Edition User Guide ([10000107224](#)) or QX200 Droplet Reader and QX Manager Software Standard Edition User Guide ([10000107223](#)) or the QX ONE Software Regulatory Edition User Guide ([100000116656](#)) for a detailed explanation of merged wells.

Section 9 Result Interpretation

The VeriCheck HEK293 Res DNA Size Kit measures two targets in the HEK293 genome. The FAM channel detects fragments above and below 200 bp and the HEX channel detects fragments ≥ 200 bp. In addition, the relative efficiencies of the two fragments are compared to estimate the percentage of all fragments larger than 200 bp as well as the median fragment length.

Qualitative Assessment

In unfragmented DNA, the two amplicons will have approximately equal size signal in the FAM and HEX channels. As DNA fragments, the ≥ 200 bp cluster in the HEX channel will shrink more rapidly than the FAM channel. Samples with no signal in the HEX channel have no detectable fragments ≥ 200 bp. Samples with no positive droplets have no detectable HEK293 DNA. See Figure 6.

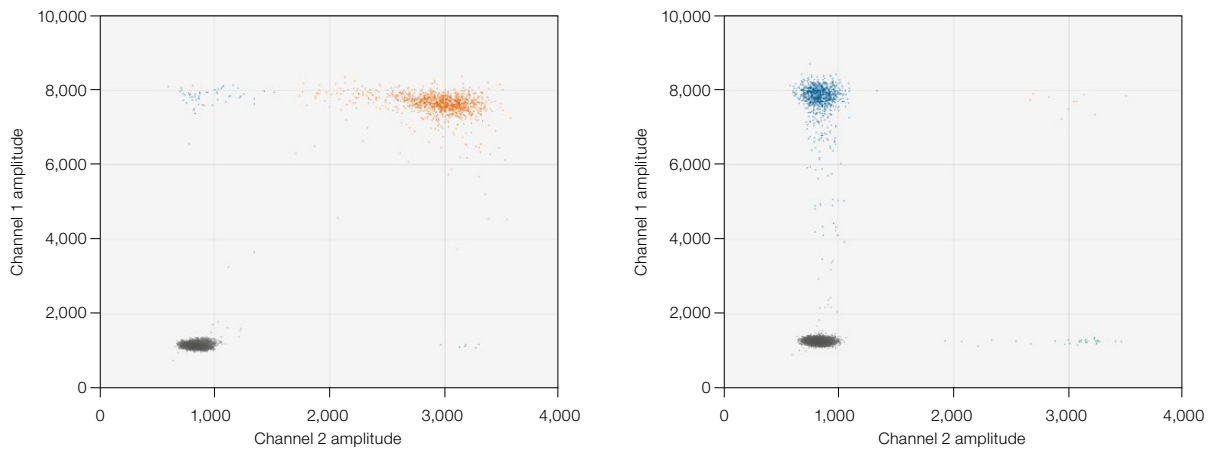


Fig. 6. 2D plot comparison of unfragmented DNA (left) versus fragmented DNA (right).

Quantitative Assessment

An additional workbook is used to quantitatively measure the level of fragmentation. To use it:

1. To export data, navigate to the Data Table tab, click the Table Menu icon on the far right, and click **Export to CSV** from the dropdown menu (Figure 7). Ensure replicate wells are merged before exporting data.

Well	Molecular Weight(ppμL)	Status	Status Reason	Experiment	SampleType	TargetType	Supernix	DyeName(s)	Copies/20μLWell	TotalConfMax	TotalConfMin	PoissonConfMa	Copy
M01	40			RDQ	PosCtrl	Unknown	ddPCR Sup...	FAM		140	124		
M01	47.3			RDQ	PosCtrl	Unknown	ddPCR Sup...	HEX		133	127		
M02	0			RDQ	Unknown	Unknown	ddPCR Sup...	FAM	0.0543	0	0.0543		
M02	0			RDQ	Unknown	Unknown	ddPCR Sup...	HEX	0.0543	0	0.0543		
M03	48.9			RDQ	Unknown	Unknown	ddPCR Sup...	FAM		138	131		138 131
M03	47.3			RDQ	Unknown	Unknown	ddPCR Sup...	HEX		133	127		133 127
M04	0			RDQ	Unknown	Unknown	ddPCR Sup...	FAM	0.0569	0	0.0569		0
M04	0			RDQ	Unknown	Unknown	ddPCR Sup...	HEX	0.0569	0	0.0569		0
M05	46.7			RDQ	Unknown	Unknown	ddPCR Sup...	FAM		135	123		132 125
M05	45.7			RDQ	Unknown	Unknown	ddPCR Sup...	HEX		140	114		129 122
M06	0			RDQ	Unknown	Unknown	ddPCR Sup...	FAM	0.0563	0	0.0563		0
M06	0			RDQ	Unknown	Unknown	ddPCR Sup...	HEX	0.0563	0	0.0563		0
M07	51.8			RDQ	Unknown	Unknown	ddPCR Sup...	FAM		150	135		146 139
M07	50.3			RDQ	Unknown	Unknown	ddPCR Sup...	HEX		147	130		142 135
M08	0			RDQ	Unknown	Unknown	ddPCR Sup...	FAM	0.0572	0	0.0572		0
M08	0			RDQ	Unknown	Unknown	ddPCR Sup...	HEX	0.0572	0	0.0572		0
M09	45.4			RDQ	Unknown	Unknown	ddPCR Sup...	FAM		132	117		128 122
M09	44.6			RDQ	Unknown	Unknown	ddPCR Sup...	HEX		126	120		126 120
M10	0			RDQ	Unknown	Unknown	ddPCR Sup...	FAM	0.0576	0	0.0576		0
M10	0			RDQ	Unknown	Unknown	ddPCR Sup...	HEX	0.0576	0	0.0576		0
M11	48.3			RDQ	Unknown	Unknown	ddPCR Sup...	FAM		137	128		136 130
M11	47.2			RDQ	Unknown	Unknown	ddPCR Sup...	HEX		133	127		133 127

Fig. 7. Screenshot of data export example.

- Open the workbook and use green **Insert QX data** button.
- Select the exported data to analyze and click **Open**.

In the result sheet you will see the following metrics, explained in Table 8.

Table 8. Explanation of metrics.

Metric	Explanation
Median size estimate, bp	Median estimated fragment size in the sample based on amplified targets
Median estimate confidence interval (CI)	95% CI for the median estimated fragment size
Percentage >200 bp	Percentage of the input sample with fragments ≥200 bp in size, estimated based on ddPCR data
Percentage >200 bp CI	Confidence interval for the estimation of the percentage >200 at 95%*
Mass, pg/μl	Calculated molecular weight from QX Software, based on amplified targets

* The values represented by the estimates of median size and percentage >200 are based on comparison of ddPCR data to quantitative gel-based sizing methods.

The confidence interval (CI) refers to error based on the detectable amplicons in the HEX and FAM channels. At low DNA concentrations, confidence intervals for these estimates will be larger. Based on the values presented in these columns, it is possible that the estimated median size is <200 bp but the CI shows the upper bound is >200 bp. In this case discretion is required to determine whether further fragmentation is required. Figure 8 shows results of an example data analysis.

VeriCheck ddPCR Residual DNA Sizing Results: HEK293											
Well	Sample Description 1	Sample Description 2	Sample Description 3	Sample Description 4	Median Size < 200 bp	Median Size Estimate, bp	Median CI, bp	Percent > 200 bp	Percent > 200 bp CI	pg / μl FAM Channel	pg / μl HEX Channel
M01	S1				No	>2000	>2000 - >2000	100%	100 - 100%	599	497
M02	S1				No	1980	784 - >2000	98.62%	89.43 - 100%	5.03	4.38
M03	S2				No	1165	1075 - 1267	93.1%	92.34 - 94.14%	570	386
M04	S2				No	1197	623 - >2000	93.62%	85.66 - 100%	6.48	4.43
M05	S3				No	880	821 - 946	89.67%	88.74 - 90.62%	605	374
M06	S3				No	577	366 - 1061	83.01%	74.57 - 93.88%	6.67	3.56
M07	S4				No	605	571 - 648	84.01%	83.01 - 85.04%	564	300
M08	S4				No	617	369 - 1285	84.42%	75.21 - 96.32%	5.41	2.9
M09	S5				No	235	229 - 244	57.67%	56.28 - 59.22%	459	122
M10	S5				No	228	177 - 311	56.09%	44.41 - 72.88%	5.38	1.37
M11	S6				No	212	207 - 218	52.4%	50.93 - 54.01%	448	104
M12	S6				Yes	191	147 - 259	46.09%	31.43 - 70.4%	3.85	0.756

Fig. 8: Example results from data analysis workbook.

Section 10 Appendices

Appendix A. Performance Specifications

Limit of Blank/Limit of Detection

The kit was tested with a negative control (water). Zero positive droplets were observed in $\geq 95\%$ of tests, supporting a limit of blank (LOB) of zero. A low concentration sample was then serially diluted in 1:1 steps. The lowest concentration sample detected above the LOB in $\geq 95\%$ of tests was at approximately 2 pg/well, the limit of detection (LOD). These data are summarized in Table 9.

Table 9. Limit of blank/limit of detection data summary.

Input, pg/well	Original Sample, pg/ μ l	Target	n	n Detected Above LOB	Detection Rate, %	Mean Measured Copies per Test
12.7	1.59	FAM	24	24	100	20.6
		HEX	24	24	100	18.7
6.8	0.85	FAM	24	24	100	10.2
		HEX	24	24	100	10.2
3.6	0.45	FAM	24	23	96	5.3
		HEX	24	24	100	5.1
1.6	0.2	FAM	24	24	100	2.1
		HEX	24	23	96	2.6
0.8	0.1	FAM	24	15	63	1.0
		HEX	24	17	71	1.3

Specificity

The kit does not detect DNA from *E. coli*, Chinese hamster ovary (CHO), normal human, or Vero cells with $>99.99\%$ specificity.

Reproducibility

Six samples were tested across multiple instruments, operators, days, and kit lots. The coefficient of variation (CV) for all tests fell below 10% for all samples except AAV, in which concentration of the HEX amplicon was low (2.58 cp/ μ l). Reproducibility data are summarized in Table 10.

Table 10. Reproducibility data summary.

Sample	n	Target	Mean Concentration, cp/ μ l	CV, %
HEK293 Sizing Positive Control – high	24	FAM	774.27	4.30
		HEX	782.71	4.15
HEK293 Sizing Positive Control – medium	24	FAM	78.30	4.39
		HEX	79.08	3.77
HEK293 Sizing Positive Control – low	24	FAM	7.99	7.08
		HEX	8.16	7.29
HEK293 Sizing Positive Control with PBS and HSA	24	FAM	79.38	4.88
		HEX	80.41	4.85
AAV	24	FAM	90.25	4.58
		HEX	2.58	12.81
Fragmented HEK293 DNA	24	FAM	306.74	4.39
		HEX	137.04	4.97

AAV, adeno-associated virus; CV, coefficient of variation; HSA, human serum albumin; PBS, phosphate buffered saline.

Dynamic Range

The kit is linear throughout the dynamic range from 2 pg/well to 300 ng/well. Fragment size can be accurately determined in samples ranging from 25 pg/well to 300 ng/well (Figure 9).

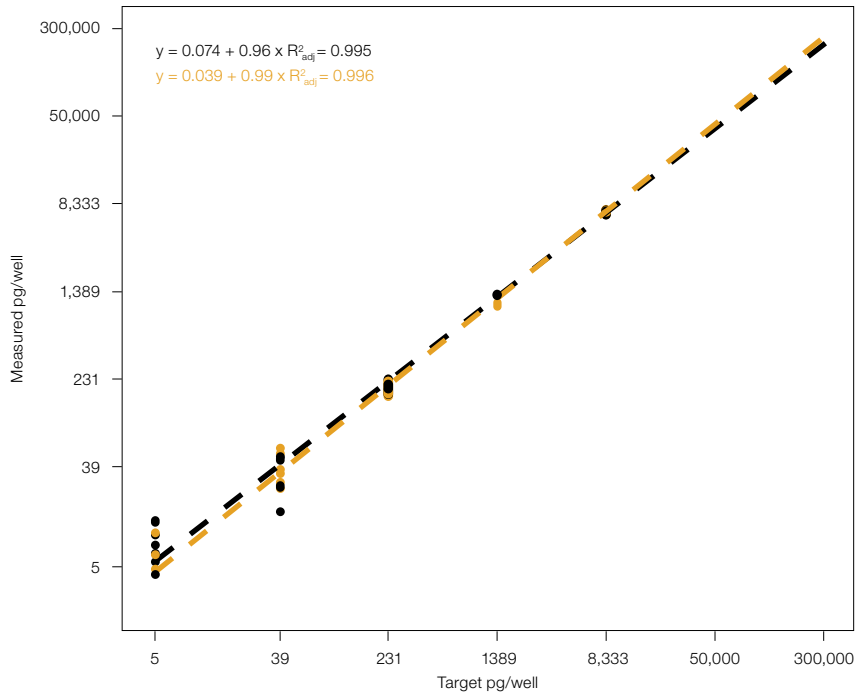


Fig. 9. Linearity. FAM channel ●; HEX channel ●.

Appendix B. Inhibitors and Extraction-Free Workflow

Samples containing ddPCR inhibitors may require dilution to meet the concentrations listed in Tables 11–13. Using samples with greater than the allowable concentration will impact results. Given an input of 8 μ l per 20 μ l reaction, the concentration in the sample is 2.5x the concentrations listed in Tables 11–13. Water is recommended for dilution, but other buffers containing low levels of Tris-EDTA, Poly(A) (Sigma-Aldrich Corporation), or Pluronic F-68 Non-Ionic Surfactant (Thermo Fisher Scientific Inc.) are also compatible. Samples containing multiple sources of inhibitors may require further dilution. Inhibitors were tested individually and not in combination with other inhibition agents. The concentrations listed in Tables 11–13 are from extraction-free samples.

Table 11. Chemical inhibitors.

Inhibitor	Maximum Allowable Concentration in Well
SDS	0.01%
Tween 20	0.01%
NP-40	0.01%
Pluronic F-68	0.625%
PEI	108 μ g/ μ l
Triton X-100	0.01%
EDTA	1.25 mM
DMSO	2.5%

DMSO, dimethyl sulfoxide; EDTA, ethylenediaminetetraacetic acid; PEI, polyethyleneimine; SDS, sodium dodecyl sulfate.

Table 12. Workflow/in-process inhibitors.

Inhibitor	Maximum Allowable Concentration in Well
EMEM Complete Media	20%
NaCl*	26 mM
Lysate (0.2% SDS)	1.25%
Histidine, sucrose, and Tween 20	5.6 mM, 93.1 mM, and 0.01%

*NaCl was formulated with a Tris buffer and not tested as an individual component.

Table 13. Sample inhibitors.

Inhibitor	Maximum Allowable Concentration in Well
AAV	>1.25e12 GC/ml
PBS + 2.5% HSA	15%

AAV, adeno-associated virus; GC, gene copies; HSA, human serum albumin; PBS, phosphate buffered saline.

Appendix C. Understanding Conversion Factor

The conversion factor (CF) provides a link between observed copies of DNA and the total mass of DNA. The conversion factor is calculated by comparing the copy number of the HEK293-specific region with a stable genomic reference. For HEK293A and HEK293T, this should be set to 0.55 cp/pg. HEK293 is a pseudotriploid cell line with a total genome mass of approximately 9.6 pg. For a stable target on a trisomic chromosome,

$$CF = \frac{3 \text{ copies}}{\text{genome}} \times \frac{\text{genome}}{9.6 \text{ pg}} = 0.31$$

To calculate the conversion factor for another HEK293 line, test an identical sample with the VeriCheck HEK293 Res DNA Size Kit, an *RPP30* reference assay (dHsaCP2500350), a *CRCP* reference assay (dHsaCP2506271), and an *MRPS5* reference assay (dHsaCP2506282), with nine wells of each assay and approximately 10 ng of DNA per well. Samples should have only clean HEK293 DNA; background human DNA will interfere with measurement.

All assays should be prepared and thermal cycled according to the package insert and digested in droplets by adding 0.5 µl of HindIII restriction enzyme directly to the master mix.

Calculate the average concentration for each assay and the CF using the equation below. If one of *RPP30*, *CRCP*, or *MRP25* is 30% lower or 25% higher than the other two, that assay should be discarded, as it indicates karyotypic drift. In this case, calculate the average of the remaining two assays.

$$CF = \frac{3 \times 0.31 \times \text{concentration}_{(\text{HEKsizing})}}{(\text{concentration}_{RPP30} + \text{concentration}_{CRCP} + \text{concentration}_{MRPS5})}$$

Section 11

Troubleshooting

This section lists some common failure modes with their phenotypes, descriptions, and suggested resolution. For a complete list of failure modes, refer to the Droplet Digital PCR Applications Guide ([bulletin 6407](#)) and the instruction manual of the instrument.

Too Much Unfragmented DNA

Problem: No distinct clusters in the 2D plot and droplet counts <10,000 (Figure 10). Wells have a high background of unfragmented DNA.

Resolution: Add 0.5 μ l HindIII restriction enzyme (10,000 units/ml) to each 20 μ l ddPCR reaction to digest DNA within droplets. In this case, add only 7.5 μ l sample to ensure a total volume of 20 μ l per reaction.

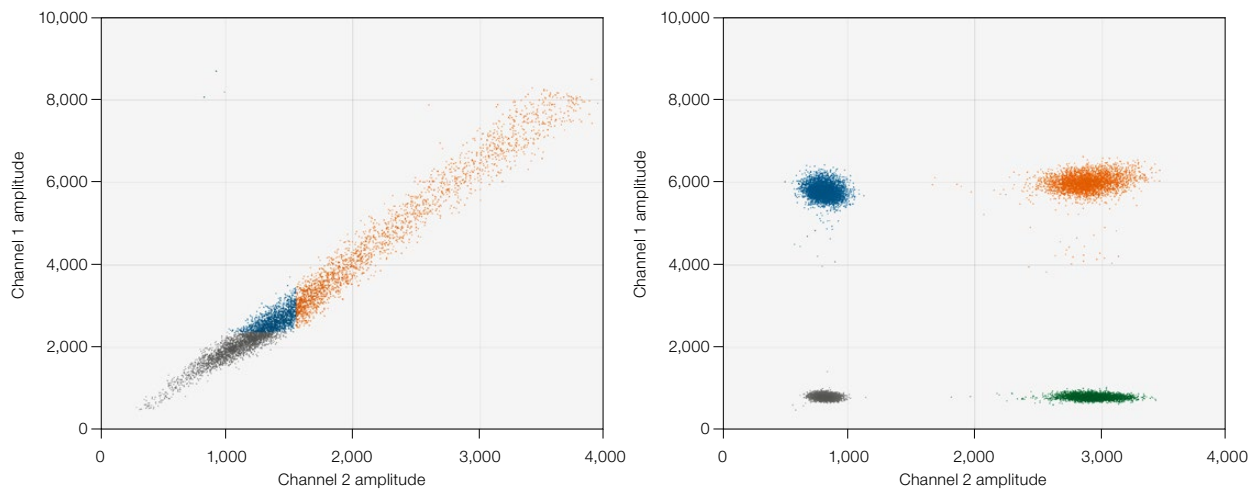


Fig. 10. 2D plots for wells containing 1 μ g of normal human DNA and 10 ng HEK293 DNA. **Left**, problem well where DNA is not digested versus **right**, normal well where restriction enzyme was added to master mix to digest background DNA.

No Negative Droplets

Problem: Only high-amplitude clusters are present in sample well. The sample is too concentrated (Figure 11) and is outside of the ddPCR dynamic range.

Resolution: Exclude well from analysis, dilute sample, and test again.

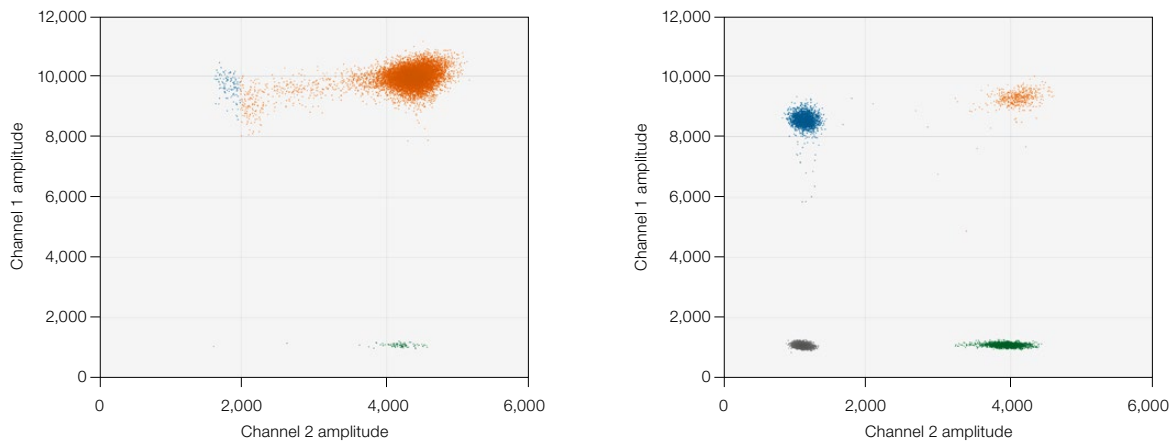


Fig. 11. No negative droplets. Left, problem well versus right, normal well.

No or Low Total Droplet Counts in Sample Wells Only

Problem: Multiple replicate sample wells have a droplet count <10,000 (Figure 12). If control wells have droplet counts >10,000, the issue is likely caused by inhibitors carried over from the starting matrix or extraction.

Resolution: Exclude sample wells with low droplet counts from analysis. If all replicate sample wells have low total droplet counts, dilute the sample to decrease inhibitor concentration and repeat Droplet Digital PCR. Alternatively, extract the sample to eliminate inhibitors.

No or Low Total Droplet Counts in Both Sample and Control Wells

Problem: Sample and control wells have droplet counts <10,000 (Figure 12). If control wells have low droplet counts, inhibition is not the likely cause of the low counts.

Resolution: Exclude wells with low total droplet counts from analysis. Repeat Droplet Digital PCR and wait 30 minutes or more at 4°C after thermal cycling completes to increase droplet numbers. If still a problem, call Technical Support.

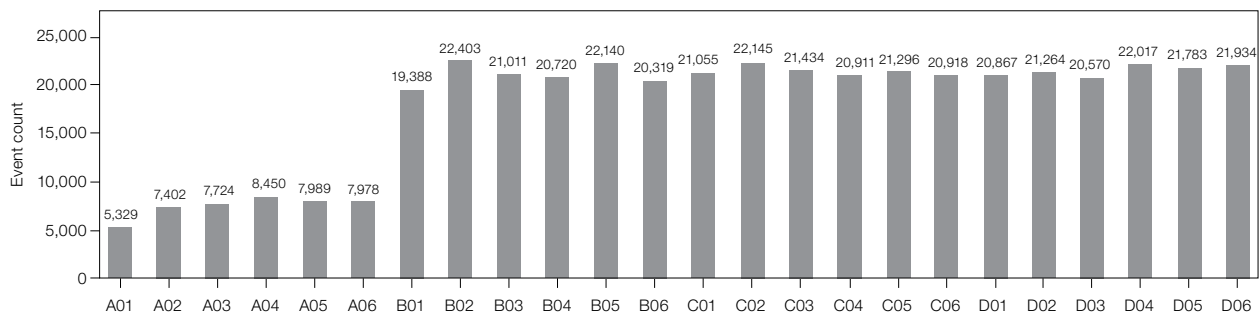


Fig. 12. Event count graph in which wells have a total droplet count <10,000.

PCR Inhibition

Problem: Separation between clusters decreases and rain increases (Figure 13).

Resolution: If concentration cannot be calculated, exclude the well from analysis. Dilute the sample to decrease inhibitor concentration and repeat Droplet Digital PCR. Alternatively, extract the DNA sample to eliminate PCR inhibitors.

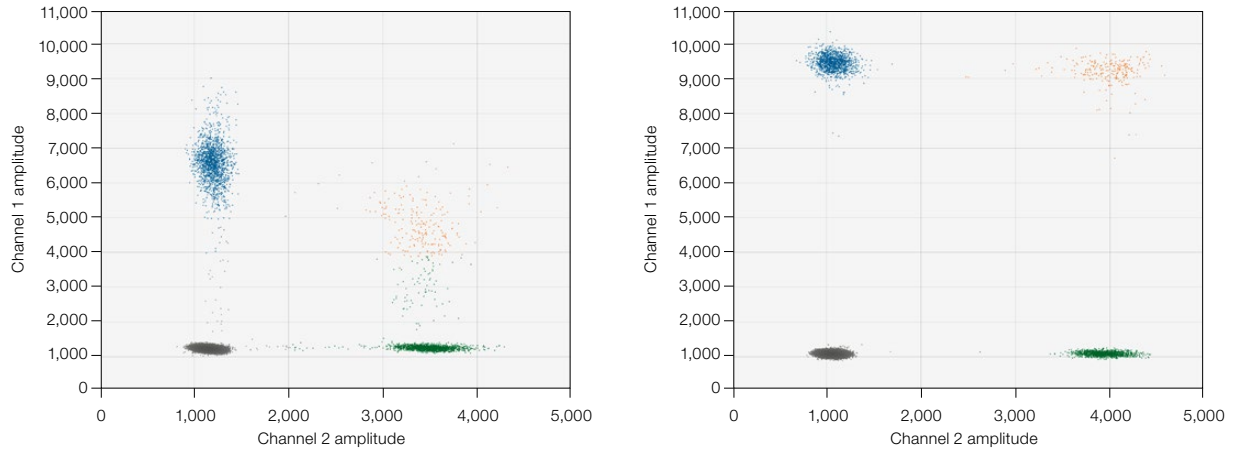


Fig. 13. PCR inhibition. Left, problem well versus right, normal well.

Mirroring

Problem: Droplets exhibit two distinct sizes (Figure 14). This indicates a potential consumable failure or particulates from samples, environment, tips, or reagents.

Resolution: Exclude the well from analysis and repeat Droplet Digital PCR, preferably with a different lot of droplet generation consumables.

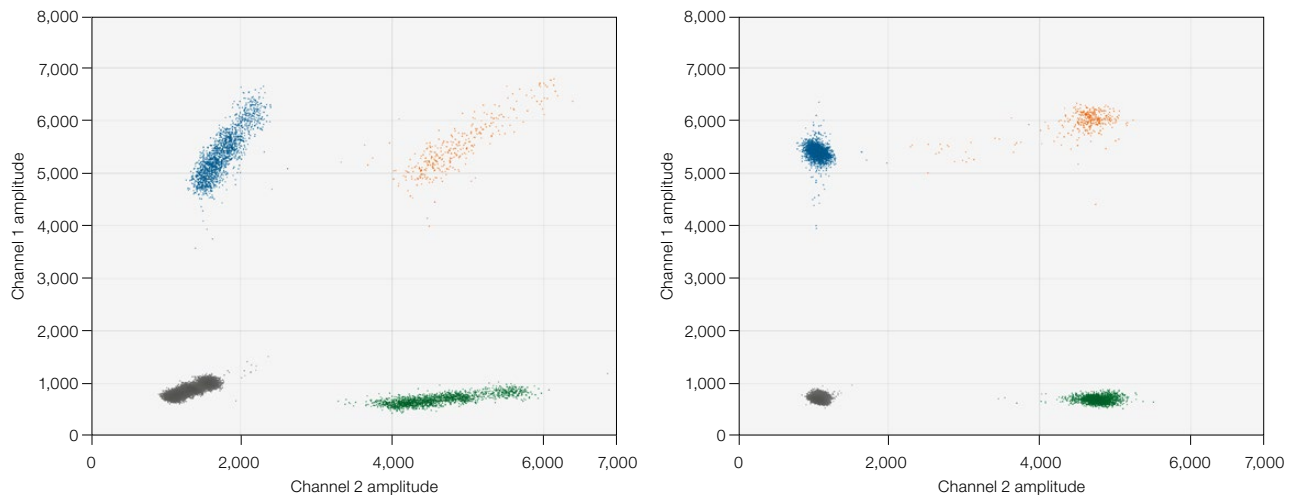


Fig. 14. Mirroring. Left, problem well versus right, normal well.

Droplet Shredding

Problem: Shredded droplets appear on the diagonal through the negative droplet cluster (Figure 15).

Resolution: Exclude well from analysis or repeat Droplet Digital PCR.

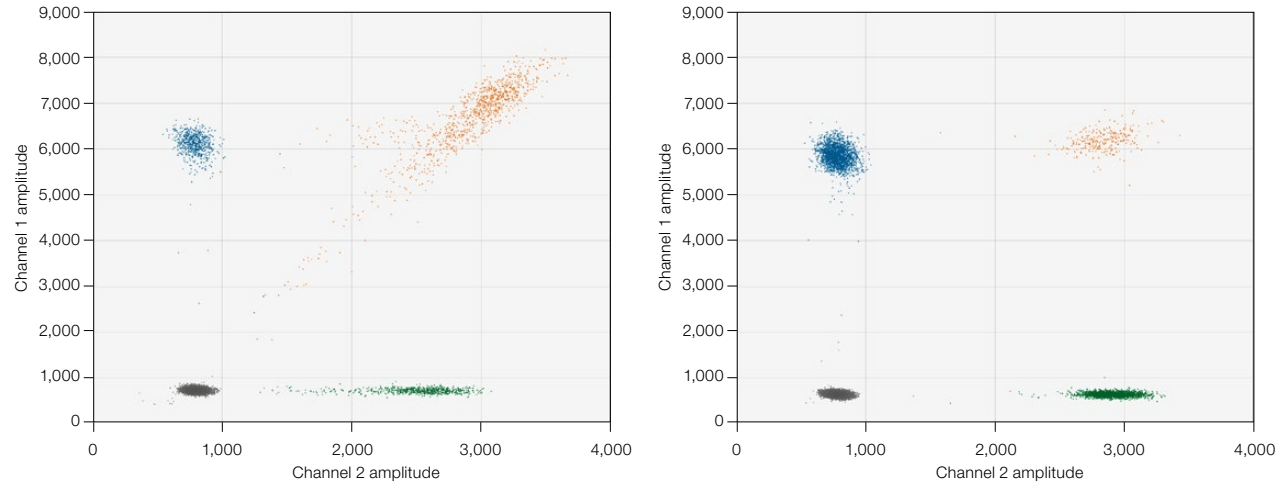


Fig. 15. Droplet shredding. Left, problem well versus right, normal well.

Incorrect Autothresholding

Problem: Droplets are incorrectly thresholded using autothresholding (Figure 16).

Resolution: Ensure that wells were thresholded using the positive control–based autothresholding with tilt. If this does not resolve the issue, follow instructions for manual thresholding. If necessary, refer to the Droplet Digital PCR Applications Guide ([bulletin 6407](#)) for best practices setting manual thresholds.

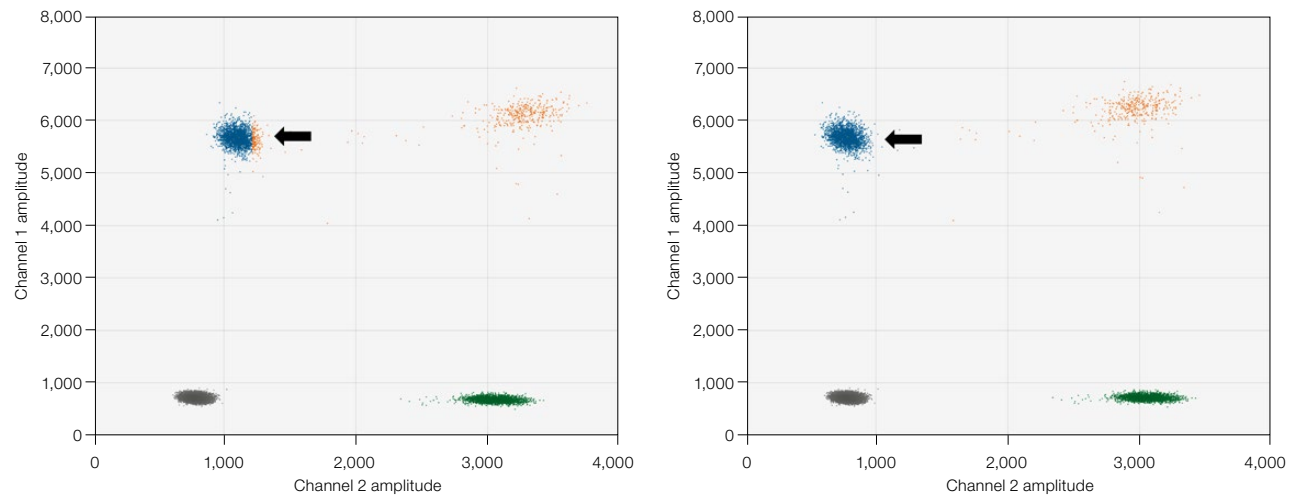


Fig. 16. Incorrect autothresholding, indicated by arrows. Left, problem well versus right, normal well with manually adjusted thresholds.

Visit [bio-rad.com/ddPCR-VeriCheck-HEK-Size](https://www.bio-rad.com/ddPCR-VeriCheck-HEK-Size) for more information.

BIO-RAD, AUTODG, DDPCR, DROPLET DIGITAL, and QX ONE are trademarks of Bio-Rad Laboratories, Inc. in certain jurisdictions. All trademarks used herein are the property of their respective owner. © 2022 Bio-Rad Laboratories, Inc.

Purchase of Digital PCR and/or Single-Cell NGS Sample Preparation products (the "Products") from Bio-Rad Laboratories is subject to Bio-Rad Laboratories, Inc. Standard Terms and Conditions of Sale, which can be accessed at <https://www.bio-rad.com/en-us/terms-conditions>. Unless we expressly state otherwise in additional Terms and Conditions, no rights are granted for you to distribute or resell the Products. Unless we expressly state otherwise in additional Terms and Conditions, no rights are granted for the development or commercialization of diagnostic assays for use with the Products without a license from Bio-Rad. The Products and/or their use are covered by U.S. and foreign patents and/or pending patent applications owned by or under license to Bio-Rad Laboratories, Inc. See <https://www.bio-rad.com/en-us/trademarks>.



**Bio-Rad
Laboratories, Inc.**

Life Science
Group

Website [bio-rad.com](https://www.bio-rad.com) **USA** 1 800 424 6723 **Australia** 61 2 9914 2800 **Austria** 00 800 00 24 67 23 **Belgium** 00 800 00 24 67 23 **Brazil** 4003 0399 **Canada** 1 905 364 3435 **China** 86 21 6169 8500 **Czech Republic** 00 800 00 24 67 23 **Denmark** 00 800 00 24 67 23 **Finland** 00 800 00 24 67 23 **France** 00 800 00 24 67 23 **Germany** 00 800 00 24 67 23 **Hong Kong** 852 2789 3300 **Hungary** 00 800 00 24 67 23 **India** 91 124 4029300 **Israel** 0 3 9636050 **Italy** 00 800 00 24 67 23 **Japan** 81 3 6361 7000 **Korea** 82 2 3473 4460 **Luxembourg** 00 800 00 24 67 23 **Mexico** 52 555 488 7670 **The Netherlands** 00 800 00 24 67 23 **New Zealand** 64 9 415 2280 **Norway** 00 800 00 24 67 23 **Poland** 00 800 00 24 67 23 **Portugal** 00 800 00 24 67 23 **Russian Federation** 00 800 00 24 67 23 **Singapore** 65 6415 3188 **South Africa** 00 800 00 24 67 23 **Spain** 00 800 00 24 67 23 **Sweden** 00 800 00 24 67 23 **Switzerland** 00 800 00 24 67 23 **Taiwan** 886 2 2578 7189 **Thailand** 66 2 651 8311 **United Arab Emirates** 36 1 459 6150 **United Kingdom** 00 800 00 24 67 23

